Thermo Scientific ARL EQUINOX 100

Transportable Benchtop X-Ray Diffractometer

The Thermo Scientific[™] ARL[™] EQUINOX 100 transportable benchtop X-Ray diffractometer provides very fast measurements compared to other diffractometers due to our unique EQUINOX curved detector that can measure all diffraction peaks simultaneously across a wide angular range. It will take you just a few minutes on most samples to get a complete analysis regardless of your resolution requirements.



Environment



Metallurgy





No scanning required. Measures the whole 2 theta range simultaneously and in real time. No external water cooling required, plugs into a standard electrical power outlet.

A rugged multi purpose benchtop style instrument designed for QA/QC, academic and routine X-ray diffraction applications.

The ARL EQUINOX 100 is the perfect instrument when space becomes a priority thanks to its water cooling system totally integrated within the instrument. With its excellent performance it is the perfect match for laboratories that are looking for a simple and powerful benchtop X-ray diffractometer. Several sample holders are available on this instrument such as a 6 positions sample holder, powder transmission mode sample holders, and a small thin film sample attachment.

Real time acquisition Easy to use No alignment needed


Very reliable Superb resolution No maintenance



Education

Specifications

X-ray source	Generator: 50 W (50kV / 1mA - Option 30kV limitation) Micro-focus X-ray tube (Cu, Co, Ag, Mo)
Detector	Curved Position Sensitive X-ray Detector, CPS 180 Curvature radius: 180mm Acquisition in real time over 110° 2Theta
Goniometric table	No moving parts: acquisition in asymmetric mode
Sample holders	Fixed, with or without sample rotation Reflection or transmission mode Zero background sample cup Controlled Atmosphere Reflection Sample Holder Gandolfi chamber Sample changer: 6 positions Thin layer attachment
Optic	Monochromatic (X-ray mirror)
Computer	Windows© 8, 7 Vista, XP
Software	Real time diffractogram display Multiple and automatic recording Peak search Deconvolution with several shapes Phase identification and quantification Degree of crystallinity determination Cell parameters, crystallite size, lattice strain Crystal structure analysis Rietveld analysis Free open database for search match Option: ICDD PDF2 or PDF4 databases
Power	Standard: 16A/ 110-230V – 50-60Hz No external water cooling
Weight	~85 kg
Dimensions	530 mm H x 520 mm D x 750 mm W
Weight	No external water cooling ~85 kg

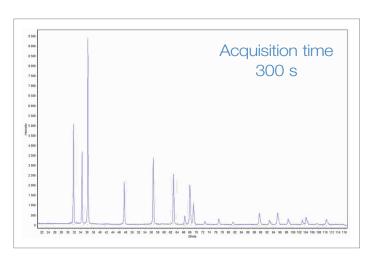

ZnO powder sample recorded on ARL EQUINOX 100 in reflection mode.



You want to perform transmission analysis? Choose either the Sample Spinner for reflection and transmission analysis on Powder samples or the Capillary sample holder for Transmission mode analysis depending on your needs.

You need to analysis several powder samples a day? Optimize your throughput with the 6 position sample changer.






You need to perform grazing incidence measurement on small samples? Try our Small Thin Layer analysis attachment.

Gandolfi chamber



Controlled atmosphere reflection sample holder

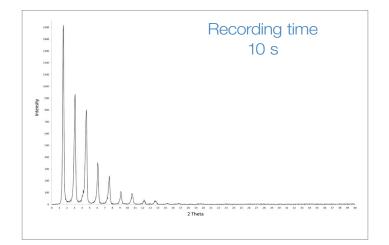

Silver Behenate powder sample recorded on ARL EQUINOX 100 in reflection mode.

www.thermoscientific.com/xrd

©2016 Thermo Fisher Scientific Inc. All rights reserved. Microsoft and Windows are registered trademarks of Microsoft Corporation. ISO is a trademark of the International Standards Organization. All other trademarks are the property of Thermo Fisher Scientific and its subsidiaries. Specifications, terms and pricing are subject to change. Not all products are available in all countries. Please consult your local sales representative for details.

Africa +43 1 333 50 34 0 Australia +61 3 9757 4300 Austria +43 810 282 206 Belgium +32 53 73 42 41 Canada +1 800 530 8447 China +86 21 6865 4588

Denmark +45 70 23 62 60 Europe-Other +43 1 333 50 34 0 Finland/Norway/Sweden +46 8 556 468 00 France +33 2 38 80 45 45 Germany +49 6103 408 1014

India +91 22 6742 9494 Italy +39 02 950 591 Japan +81 45 453 9100 Latin America +1 561 688 8700 Middle East +43 1 333 50 34 0 Netherlands +31 76 579 55 55

Thermo Fisher Scientific (Ecubiens) SARL, Switzerland is ISO certified.

New Zealand +64 9 980 6700 Russia/CIS +43 1 333 50 34 0 Spain +34 914 845 965 Switzerland +41 61 716 77 00 UK +44 1442 233555 USA +1 800 532 4752

